



FreeStyle™

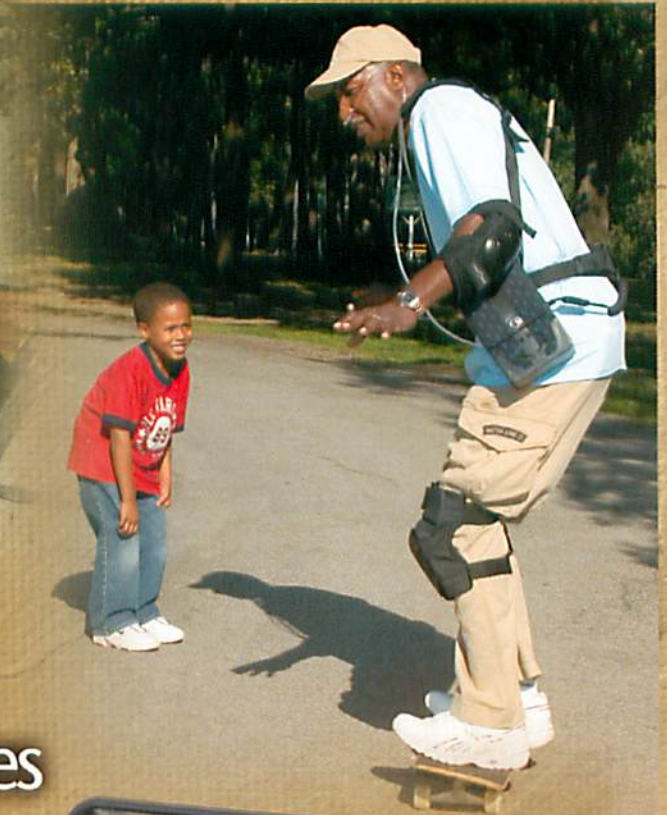

For **Oxygen Patients** with No Boundaries

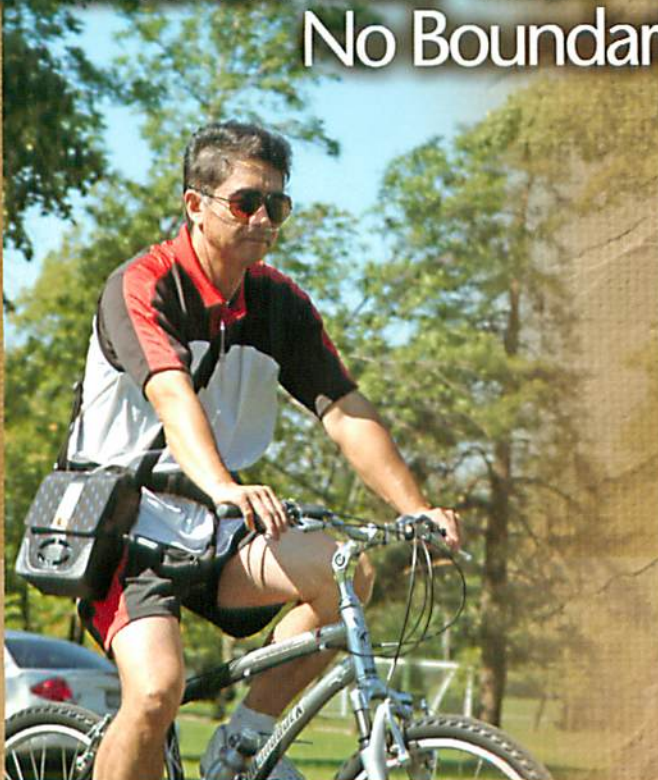

FREESTYLE™ PORTABLE OXYGEN CONCENTRATOR

AIRSEP®

FreeStyle™

Let Freedom Reign!

FreeStyle™ by AirSep is the only portable oxygen concentrator tailored specifically for those oxygen patients who know no boundaries. ...Oxygen patients who present the ultimate challenge to you and your busy staff. ...Oxygen patients who are dissatisfied with the limitations imposed by contents-based gas and liquid portable systems.

Patients like these define the essence of FreeStyle. The FreeStyle Portable Oxygen Concentrator represents a ticket to freedom for themselves and their loved ones because it is the most compact, lightweight, full-featured system to enter the supplemental oxygen arena. Bold and provocative like its name suggests, FreeStyle captures the heart of what it means to live the good life – without limitation. This portable concentrator offers the same weight advantage of a contents-based system, yet provides an unlimited supply of oxygen. Supplemental oxygen users can now participate in an all-day, outdoor music fest, a sporting event, a continuing education class, or even enjoy the added benefits of FreeStyle during extended travel.

FreeStyle Ensembles

To maximize ambulation time, FreeStyle can be used as an ensemble consisting of the portable oxygen concentrator, comfortable carrying case, and matching AirBelt. The unit provides three pulse settings equivalent to 1–3 LPM oxygen. A light above the pulse settings


illuminates to show when maintenance is required. The flexible, swivel oxygen outlet allows for comfortable positioning of the cannula. Through the transparent control panel cover, the unit also displays a battery test button, which is used to verify the status of the internal battery. Soft-touch buttons, enhanced with Braille on the unit to note the pulse settings, further hallmark this product's exceptional design and attention to detail.



Exceptional Battery Life

FreeStyle's extended battery pack fits inconspicuously in the exclusively designed AirBelt™. This battery belt offers not only fashion but function for those who won't leave home without the added assurance of additional hours on FreeStyle, which is equipped with its own battery indicators. And AirBelt is so lightweight and ergonomically correct that one can easily forget that it's being worn.


AirBelt connects effortlessly to the FreeStyle unit, providing the user up to 10 hours between recharges.

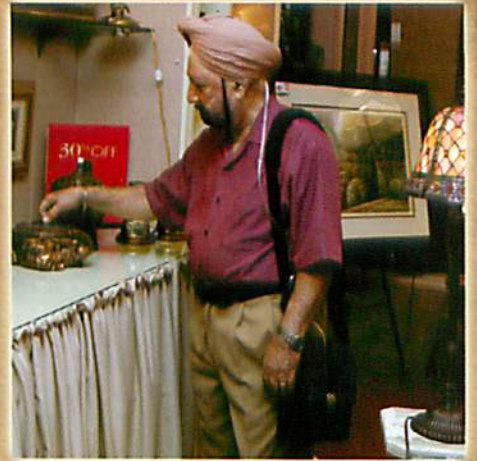

— FreeStyle — For Those with a Passion for Life

We realize that you have patients with high hopes and expectations who wish to travel more or even seek out unbridled adventure, but feel limited by their current portable oxygen systems. These are the most active of oxygen patients and those who often do not do well with limitations or restricted schedules. Struggling to make the most of limited contents based systems sometimes leads to non-compliance or running out, creating health risks ultimately affecting longevity.

Finally, with FreeStyle, you have the means to bring unlimited oxygen to their world and unleash unlimited freedom in a style you never imagined – until now.

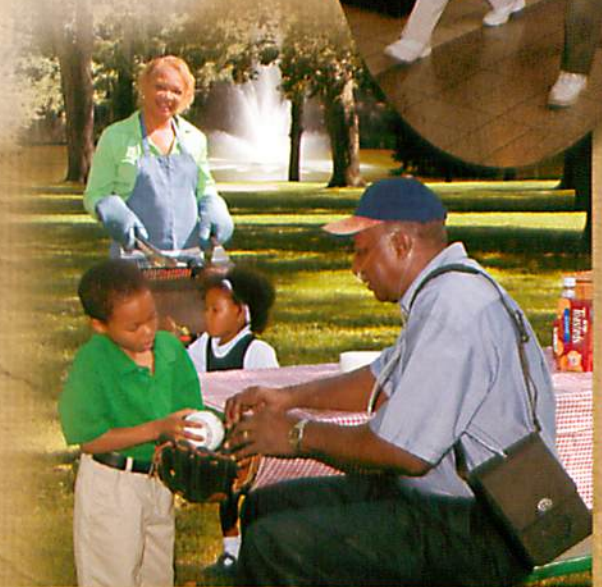

Fashion with Function

You and your patients no longer need to settle for antiquated-looking medical devices that don't perform in a manner that keeps up with the times. Modern-day patients are looking for contemporary, portable solutions that help them truly live better – and so are you. In an era of ever-increasing fuel costs, routine non-revenue deliveries tax your staff and limit your company's growth.

FreeStyle makes your patient more ambulatory, while eliminating the burden and expense of oxygen deliveries. It's a win-win situation for both you and the exceptionally active person prescribed for supplemental oxygen.

Better ways to ensure your oxygen patients' compliance is key to their continued health and longevity.

With FreeStyle at your patient's side, what was once improbable, or even impossible, is now a day-to-day reality. If your patient can do it, FreeStyle will be right there as a steady, lightweight, and long-lasting, stylish companion.



Unit and ensemble shown in designer ebony.



Unit and ensemble shown in designer beige.

SPECIFICATIONS

OXYGEN CONCENTRATION*	1–3 pulse settings equivalent to continuous flow of 1–3 LPM
DIMENSIONS:	8.6 in. H x 6.1 in. W x 3.6 in. D (21.8 cm H x 15.5 cm W x 9.1 cm D)
WEIGHT:	4.4 lb (2 kg), 1.8 lb (.8 kg) battery belt
POWER:	100–240 VAC, 50/60 Hz; 12–16 VDC
BATTERY DURATION:	Rechargeable lithium ion battery

<i>Internal Battery only</i>	3 LPM - 2 hours 2 LPM - 2.5 hours 1 LPM - 3.5 hours
------------------------------	---

<i>Internal Battery combined with external AirBelt</i>	3 LPM - 5 hours 2 LPM - 6 hours 1 LPM - 10 hours
--	--

AUDIBLE ALARMS AND VISUAL INDICATORS:	Start-up Pulse flow Battery condition Cannula disconnect System overdraw High and low pressure Maintenance required
--	---

*Based on an atmospheric pressure of 14.7 psia (101 kPa) at 70°F (21°C).

Note: Specifications above subject to change without notice.

INTERNATIONAL APPROVALS

AirSep Corporation is ISO 9001:2000-certified.

ORDERING INFORMATION

Distributed by: **intermedical**

Unit 6, Mill Hall Business Estate
Mill Hall, Aylesford, Kent
ME20 7JZ United Kingdom
T 00 44 (0)1732 522444
F 00 44 (0)1732 872883

E-mail admin@intermedical.co.uk
Website www.intermedical.co.uk



ISO 9001-2000
Certificate No.
FM 508699



ISO 13485-2003
Certificate No.
MD 508698

AirSep Corporation

Medical Division

401 Creekside Drive • Buffalo, New York 14228-2085 USA • USA/Canada Toll-Free: 800.874.0202
Tel: 716.691.0202 • Fax: 716.691.4141 • E-mail: mpd@airsep.com • Internet: www.airsep.com